


# Profiling Analysis of Polysulfide Silane Coupling Agent

Technical Report vol.42

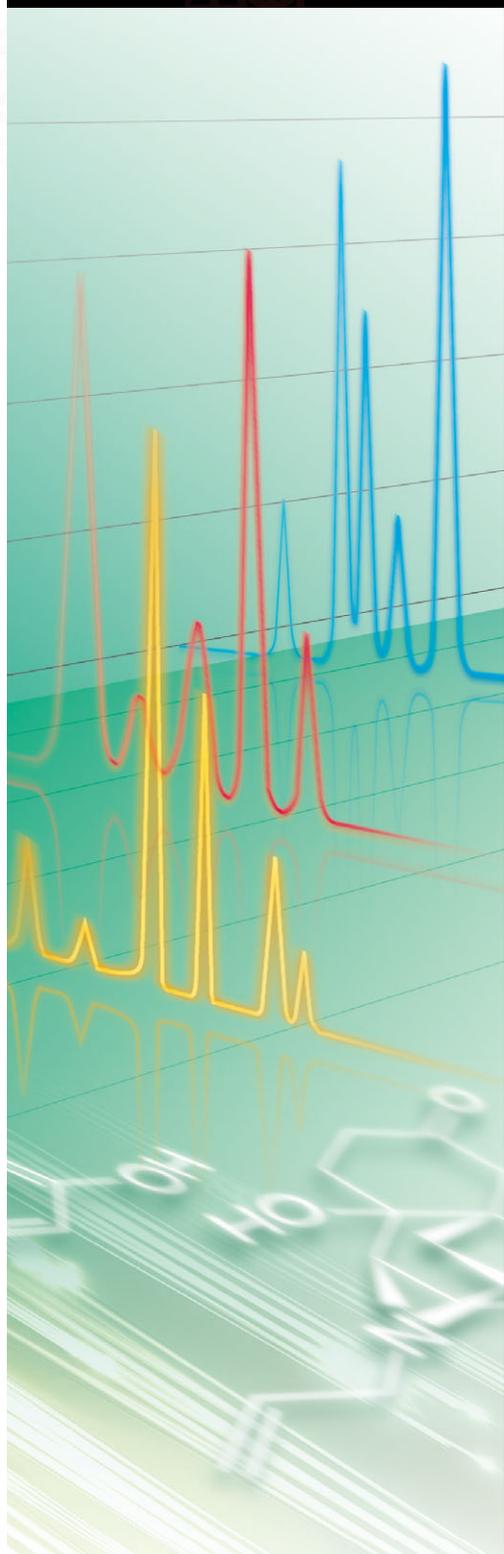

## 1. Introduction

Recently, in line with the oil-conservation movement, silica, as a non-petroleum resource, is increasingly used as a tire reinforcement filler. Fillers are typically mixed with rubber compounds used in tires to enhance reinforcement of the tire. Because silica has a hydrophilic silanol group on its surface, when added to the lipophilic elastomer rubber, it readily coheres to the rubber without actually bonding to the rubber molecule. Thus, the processibility and mechanical strength of the rubber compound is degraded, greatly affecting tire performance. For this reason, processing with a silane coupling agent is conducted to avoid this shortcoming.

Silane coupling agents contain hydrolyzable groups which react with inorganic materials, and organofunctional groups which react with organic materials. Therefore, they possess the functionality to form durable bonds between an inorganic material and an organic material. Polysulfide silane coupling agents are widely used in silica-mixed rubber compounds. The alkoxy groups at both ends of the molecule cover the surface of the silica, and the polysulfide group in the center forms a bond with the elastomer (Fig. 1). In this way, the silane coupling agent binds together the silica and elastomer to improve processibility and mechanical strength of silica-mixed rubber compounds.

Many kinds of polysulfide silane coupling agents are currently being marketed, and the silane coupling agents adopted at various tire manufacturers' development sites adhere to their own proprietary standards. Therefore, analyzing these helps to elucidate one of the characteristics associated with tire development, which is extremely important from the standpoint of tire inspection and improvement. Although it is relatively easy to analyze the differences between tires when the silane coupling agents used have completely different structures, analysis is very difficult in cases where different manufacturers are producing tires with the same chemical structure.

Assuming the analysis in the scenario above is part of a basic study, we conducted difference analysis of the principle components and sub-components in 4 different types of a commercially available polysulfide silane coupling agent, Bis (triethoxysilylpropyl) polysulfide, with the same basic skeleton from different manufacturers. The search for marker molecules that characterize the respective polysulfide silane coupling agents was conducted using Profiling Solution and multivariate analysis software. In addition, the structure analysis of the formulas and structures of the marker molecules was conducted using a combination of LCMS-IT-TOF (high performance liquid chromatograph / ion trap time-of-flight mass spectrometer), providing MS<sup>n</sup> analysis using accurate mass measurement, and Formula Predictor, a software program for prediction of molecular formulas.

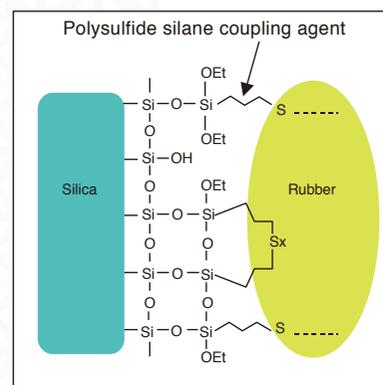


Fig.1 Silica and Elastomer Bonds Linked by Polysulfide Silane Coupling Agent

## 2. Statistical Analysis Methods

### 2-1. Multivariate Analysis

Multivariate analysis is a statistical technique for analyzing the mutual relationship between data associated with multiple variables. Using this technique helps in the identification of a specific cause associated with a particular phenomenon, and facilitates profiling analysis to determine how to group and characterize that phenomenon.

As an application example in mass spectrometry, multivariate analysis is commonly used to compare the data acquired from multiple samples, and to identify characteristic components (marker molecules) in these samples through a "marker molecule search." It is used for analyzing proteins, lipids and other types of metabolites in the body for use in the search for specific disease markers for diagnosis, as well as the search for target molecules in drug discovery. In other fields, it is an extremely useful analytical technique for finding differences and relationships that exist among multiple samples.

### 2-2. Profiling Solution

Profiling Solution is used to automatically align retention time data from the high mass accuracy data acquired using Shimadzu's LCMS-IT-TOF instrument. The results are generated in the form of a tabular matrix and can be exported (Scheme 1) to third-party multivariate analysis software (e.g. SIMCA-P; Umetrics).



Scheme 1: Multivariate Analysis Workflow Using Profiling Solution

An example of a matrix generated using Profiling Solution is shown in Fig. 2. From left to right, the columns display data files, detected ion  $m/z$  values, their retention times, and signal intensities, respectively. Besides these, multiple TIC, extracted ion chromatograms and mass spectra can also be displayed. Thus, Profiling Solution can be used to easily conduct batch processing on multiple files, making it an indispensable tool for statistical analysis.

ID	Attrib1	Filename	Ion m/z	Ion	a1	a2	a3	a4
a1	a	20-1-2kd	215.0331	2.337	11430.025	10292.818	4568.882	10292.858
a2	a	20-1-2kd	249.0620	2.654	8326.721	7.672.707	7133.551	6.737.968
a3	a	20-1-3kd	335.0518	3.064	7.845.587	8.910.588	9.003.788	9.011.696
a4	a	20-1-4kd	158.0084	2.342	6.869.401	6.700.495	6.454.390	6.231.298
			226.9507	2.704	6.933.574	7.264.561	6.178.374	5.693.545
			249.9617	3.525	6.985.289	6.678.781	6.381.980	6.654.434
			197.9092	2.308	6.751.435	5.044.388	5.016.483	6.423.326
			226.9510	3.433	4.667.503	4.963.494	4.328.902	4.112.141
			363.1182	2.433	3.720.367	3.951.726	3.942.347	3.804.764
			226.9512	3.400	3.703.856	3.729.227	3.638.912	3.425.238
			167.0214	3.060	3.373.233	3.554.626	3.088.740	3.478.484
			226.9779	7.267	3.331.592	2.932.590	2.160.700	1.856.215
			203.0533	2.383	3.227.327	3.663.151	3.949.934	3.483.838
			162.0654	2.088	3.114.621	2.690.550	1.818.344	2.411.711
			261.0145	3.021	2.969.984	1.737.425	2.432.012	2.119.054
			169.0350	3.064	2.797.278	3.076.458	3.463.383	3.301.866
			316.9484	3.403	2.515.800	2.515.421	2.500.138	2.344.126
			158.9644	2.711	2.420.310	2.485.450	2.324.742	1.925.242
			249.9630	3.038	2.282.714	3.362.336	3.018.163	1.967.308
			316.9490	2.804	2.275.459	2.102.221	1.965.951	1.812.350
			357.0344	3.038	2.253.545	2.207.268	1.227.617	1.530.511
			364.9379	3.592	2.208.700	2.318.710	3.114.781	2.103.370
			316.9480	3.433	2.201.581	2.114.736	2.338.983	2.367.982
			202.1801	7.421	2.195.347	2.578.849	2.378.481	2.256.639
			364.9375	3.417	2.133.420	2.470.783	1.860.852	1.860.568
			141.9892	2.088	2.113.332	1.517.529	1.204.700	1.325.242
			284.9389	2.729	2.063.127	2.038.494	1.746.180	1.701.583
			364.9379	2.758	2.037.531	1.841.655	1.658.718	1.442.535
			226.9780	1.142	2.014.266	1.521.921	1.364.200	1.265.302
			362.9270	2.767	1.973.263	1.903.312	1.554.129	1.535.287
			452.9268	0.382	1.974.251	1.968.142	1.694.563	1.618.780
			452.9264	3.421	1.958.565	1.645.203	1.659.423	1.855.754
			249.9611	7.648	1.925.680	1.849.601	1.567.908	1.763.378
			452.9262	2.782	1.826.983	1.440.743	1.205.961	1.158.012
			249.9615	4.167	1.897.985	1.840.286	1.521.182	1.034.570
			175.1131	3.723	1.632.448	1.346.651	1.452.342	1.524.106
			284.9383	3.423	1.610.720	1.309.412	1.655.605	1.513.548

Fig.2 Profiling Solution Window

### 3. Experimental Method

- 1) The analytical samples consisted of 4 types of a polysulfide silane coupling agent from different manufacturers which were diluted 1:10000 with methanol. In addition, the diluted solutions were mixed together at equal volumes to prepare a QA/QC sample to confirm system repeatability. Fig. 3 shows the structure of the principle component Bis (triethoxysilylpropyl) polysulfide (TESPP) of the polysulfide silane coupling agents.
- 2) LC/MS analysis was conducted using an LCMS-IT-TOF (high performance liquid chromatograph / ion trap time-of-flight mass spectrometer). A Shimadzu Shim-pack XR-ODS (75 mmL. x 2.0 mmL.D., 2.2  $\mu$ m) was used. Gradient elution analysis was conducted using a 5 mmol ammonium formate aqueous solution as mobile phase A, and acetonitrile as mobile phase B. For ionization, electrospray ionization (ESI) was used, and the positive ions were measured.
- 3) The 4 types of polysulfide silane coupling agent samples and the QA/QC sample were each analyzed 3 times. Profiling Solution was used to conduct peak picking and retention time alignment with respect to these data, and a matrix was generated. The matrix was then exported to SIMCA-P for principle component analysis (PCA), and the differences among the various polysulfide silane coupling agent compositions were expressed visually.
- 4) Many of the constituents present in the polysulfide silane coupling agents that were identified in PCA were subjected to composition prediction using formula prediction software, based on information that included accurate mass, isotope pattern, and product ion information. In addition, structure prediction based on MS<sup>n</sup> data of known substances was conducted.

### 4. Analytical Conditions

Column:	Shim-pack XR-ODS (75 mmL. x 2.0 mmL.D., 2.2 $\mu$ m)
Mobile phase A:	water containing 5 mmol ammonium formate
Mobile phase B:	acetonitrile
Time program:	70%B(0 _ 0.11 min) _ 100%B(6.67 _ 12 min) _ 70%B(12.01 _ 15 min)
Mobile phase flow rate:	0.45 mL/min
Column temperature:	40 °C
Sample injection	
Volume:	1 $\mu$ L
Ionization mode:	ESI +
Nebulizer gas flow rate:	1.5 L/min
Probe voltage:	+4.5
CDL temperature:	200 °C
BH temperature:	200 °C
Drying gas pressure:	0.1 MPa
Scan range:	<i>m/z</i> 100 - 1000

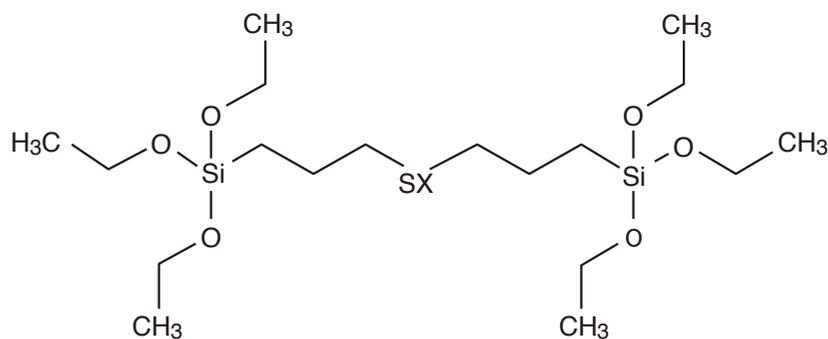


Fig.3 Structure of TESPP

## 5. Results and Discussion

### 5-1. Results of LC/MS Analysis of Polysulfide Silane Coupling Agents

Fig. 4 shows the total ion chromatograms (TICs) of the 4 types of polysulfide silane coupling agents (Sample A, B, C and D). A complete picture of the components detected using these conditions can be seen in the TIC. In addition, comparing the TICs reveals some of the differences between the samples. Multiple peaks associated with the principle components are detected, and the distribution based on the number of S was suggested. Due to the similarity of the TICs of Sample A and Sample B among the 4 sample types, the principle components in these samples were presumed to be similar overall.

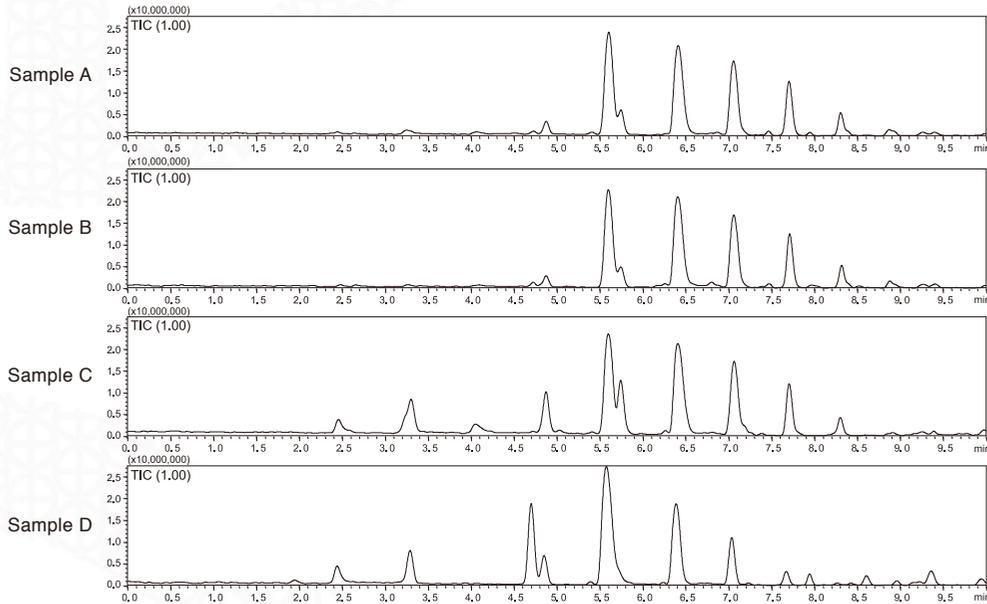


Fig.4 TICs of 4 Types of Polysulfide Silane Coupling Agents

Fig. 5 shows the mass spectra of the peaks having a retention time of 6.4 min. The ion of  $m/z$  524 was assigned to the ammonium adduct molecular ion  $C_{18}H_{42}O_6Si_2S_3$ , based on the accurate mass. This was presumed to be Bis (triethoxysilylpropyl) trisulfide (TESPT), one of the principle components of the polysulfide silane coupling agents, and the presumption was confirmed based on the observed product ion spectra obtained from the  $MS^2$  and  $MS^3$  analysis.

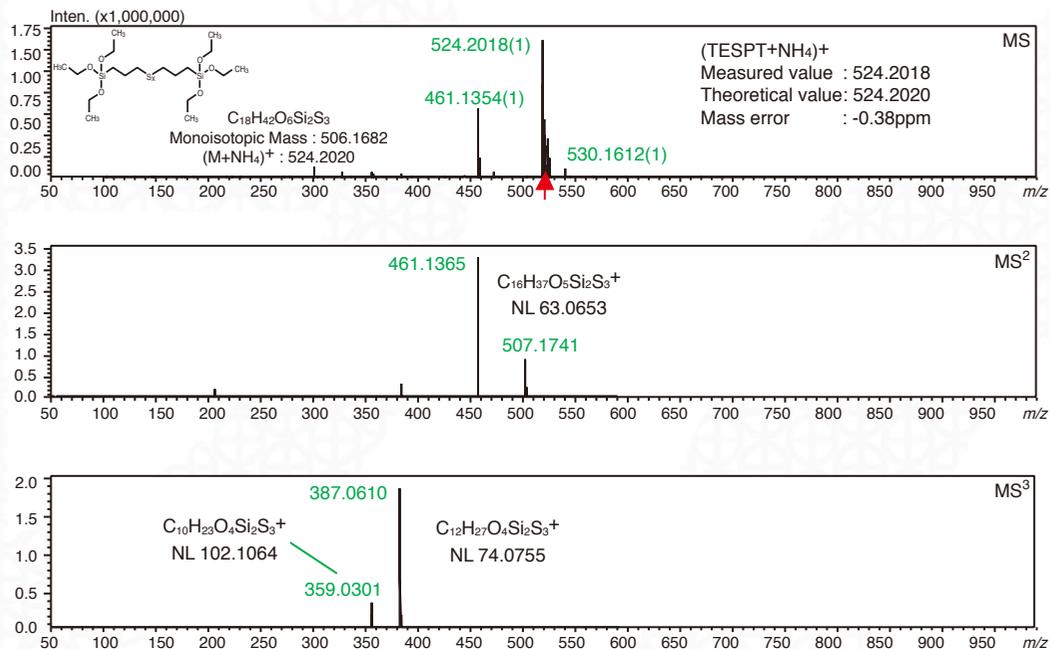


Fig.5 Mass Spectra of TESPT

Fig. 6 shows the structure of TESPT, in addition to the predicted cleavage positions of the product ions obtained from MS<sup>2</sup> and MS<sup>3</sup>. Prediction of the ion cleavage positions in the known substance based on its structure provides useful information for predicting the structures of unknown related substances.

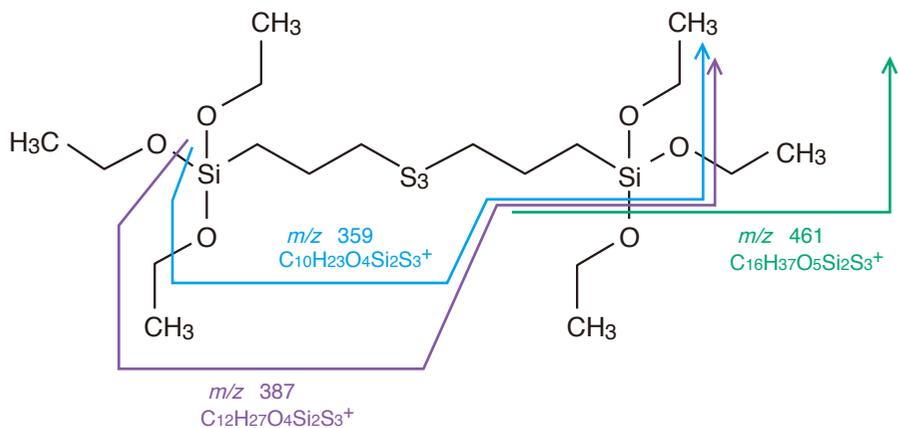


Fig.6 Predicted Cleavage Positions of TESPT

Fig. 7 shows the extracted ion chromatograms of the ammonium adduct molecular ion of TESPP (x = 1 to 8). Comparing the chromatograms makes it possible to compare the x distribution in TESPP. Similar distributions are found in Samples A, B and C, confirming that x = 2 to 8. On the other hand, the distribution in Sample D compared to that in the other 3 samples indicates a shift of x to a smaller value, indicating that x = 1. This confirmed the presence of Bis (triethoxysilylpropyl) monosulfide (TESPM) only in Sample D.

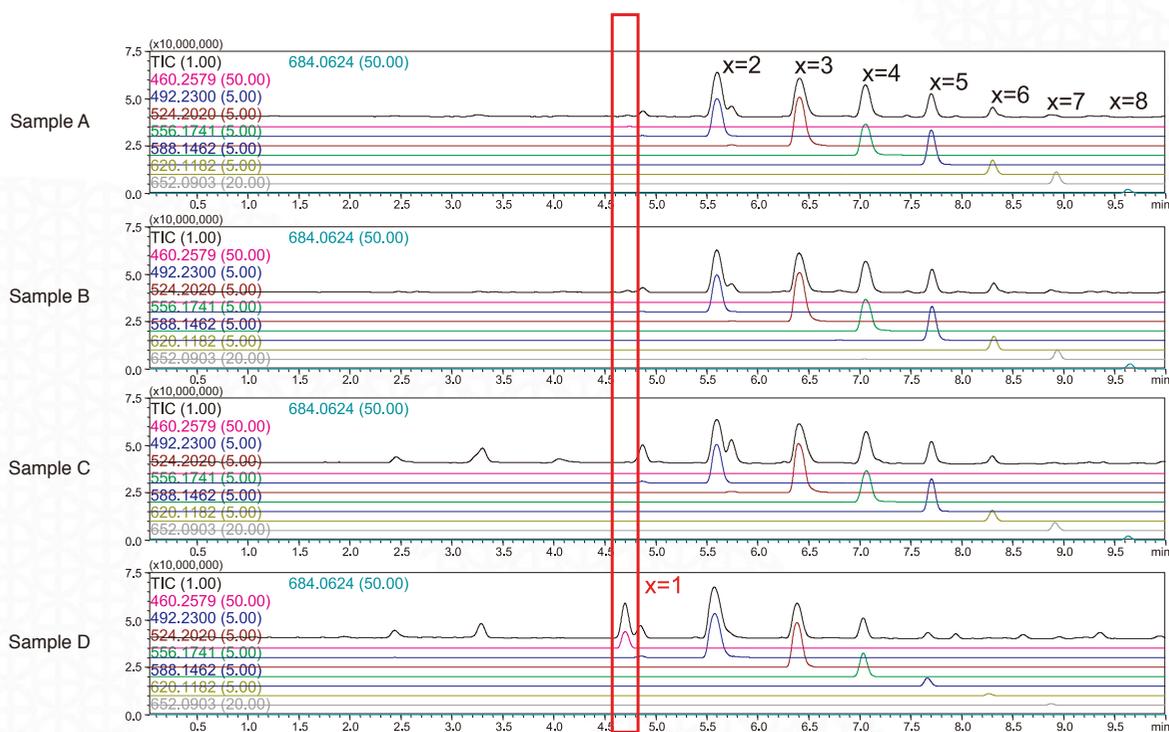


Fig.7 Extracted Ion Chromatograms of 4 Types of Polysulfide Silane Coupling Agents

The LC/MS analysis results confirmed the similarity of component compositions of Sample A and Sample B, and that the distribution of the principle component of Sample D is different from that in the other 3 samples.

## 5-2. Principle Component Analysis (PCA)

The data obtained above were used to create a matrix after picking peaks and aligning retention times using Profiling Solution (Fig. 8). In the matrix, in addition to the retention time and  $m/z$  values, the signal intensities are displayed using different colors. Also, the extracted ion chromatograms of components of interest can be displayed to check the peak shape.

The matrix generated in Profiling Solution was exported to SIMCA-P, and principle component analysis (PCA) was conducted. Fig. 9 shows the PCA results. Each point in the score plot displays 1 point of analysis data, and the closer the point positions are to one another, the higher the degree of data similarity. The high similarity of Sample A and Sample B as judged from the TIC is in agreement with the high degree of similarity indicated by their close positional distribution in the score plot (Fig. 9).

Each point in the loadings plot represents a peak in the extracted ion chromatograms, and the points that are distributed in the same direction as each of the points in the score plot are judged to have a large degree of influence on that data (indicating characteristic similarity).

In this example, the direction of distribution of Var\_355, Var\_395 and Var\_272 in the loadings plot is the direction of distribution of Sample C in the score plot, leading to the conclusion that these peaks show the characteristics of Sample C. Similarly, it was concluded that Var\_142, Var\_235, Var\_113 and Var\_284 show the characteristics of Sample D. Fig. 9 shows the list of peaks characteristic of Sample C or Sample D from the loadings plot.

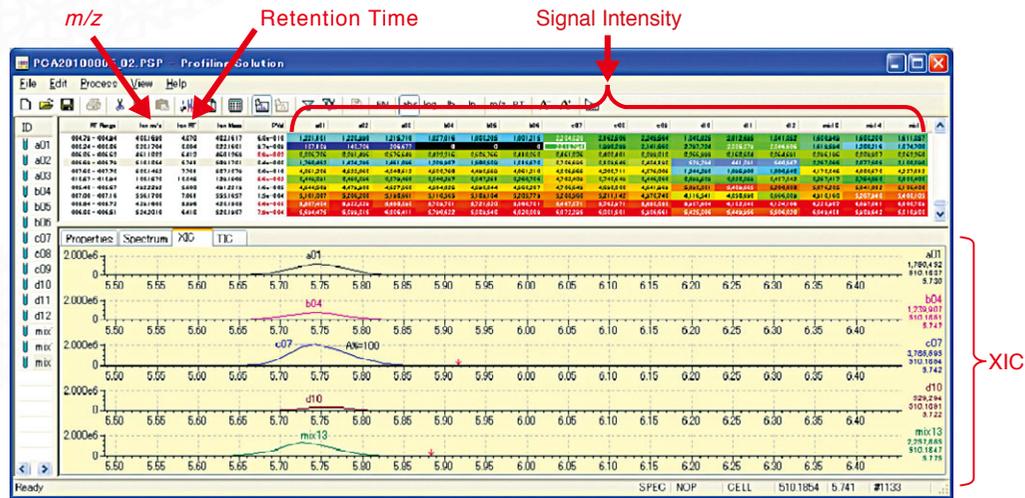


Fig.8 Matrix Generated Using Profiling Solution

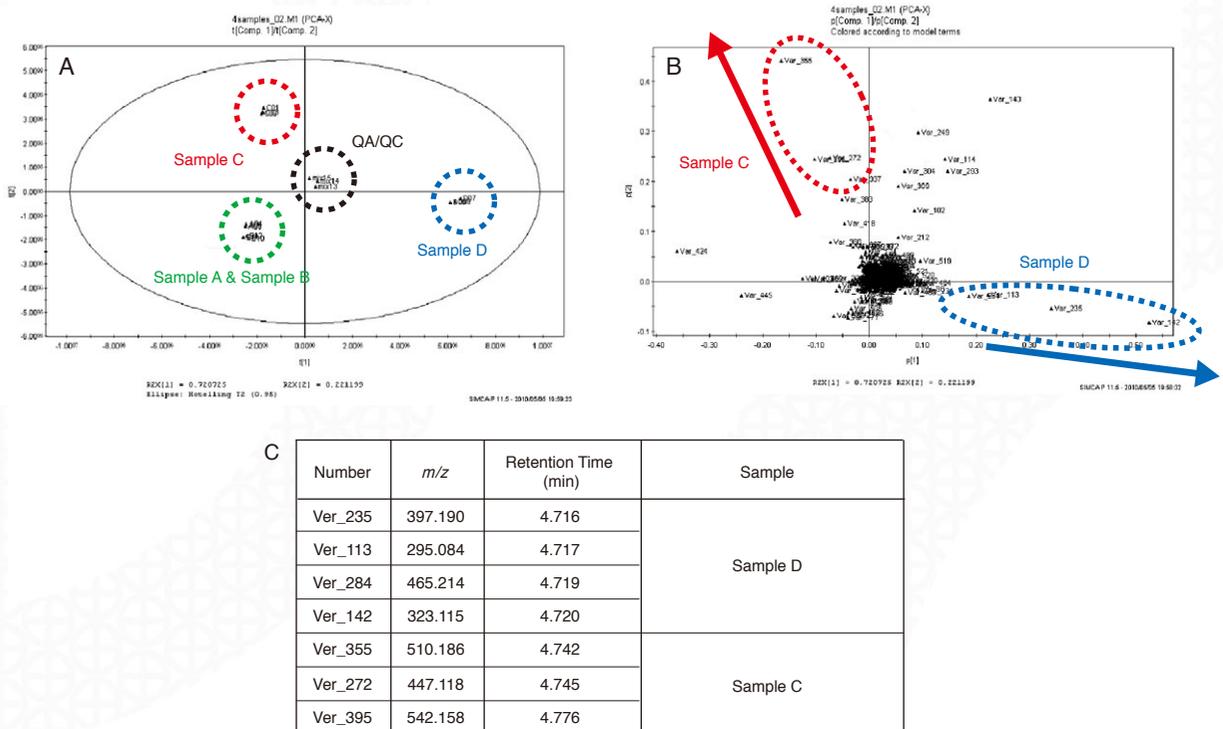


Fig.9 Principle Component Analysis (PCA) Results for 4 Polysulfide Silane Coupling Agent Samples and List of Peaks Characteristic of the Samples (A: Score Plot, B: Loadings Plot, C: List of Peaks Indicating Characteristic Association)

## 5-3. Formula and Structure Prediction

Fig. 10 shows the extracted ion chromatograms of peaks indicating characteristic associated with Sample C or Sample D. We confirmed that Component 1, derived from  $m/z$  295, 323, 397, and 465, exists in large quantity in Sample D. On the other hand, Component 2, derived from  $m/z$  510 and 447, and Component 3 derived from  $m/z$  542, were confirmed to be present in large quantity in Sample C.

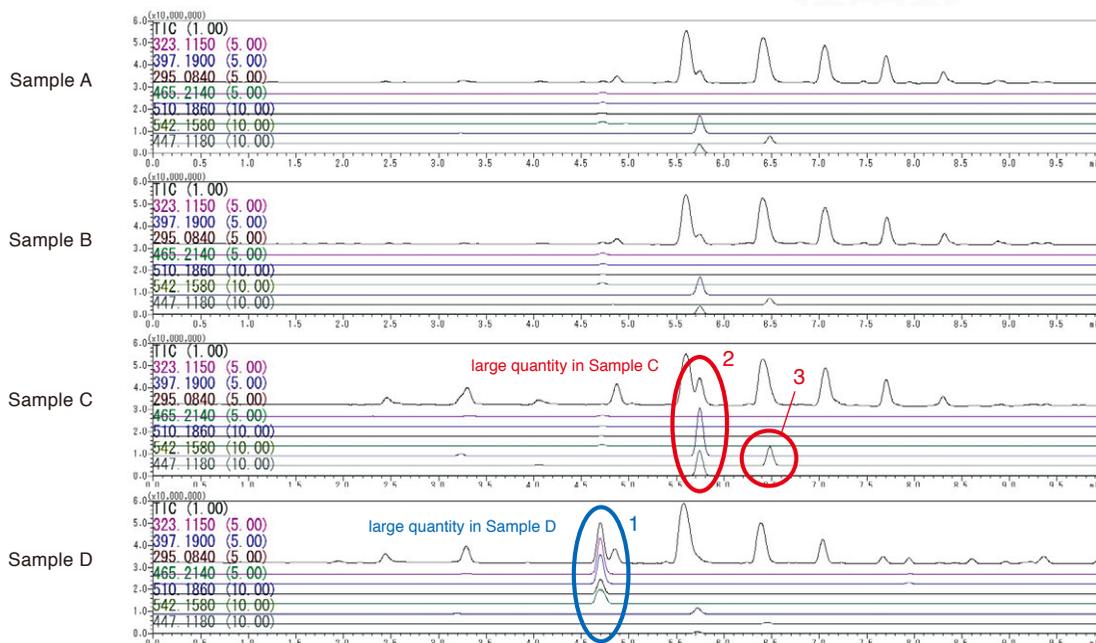


Fig.10 Extracted Ion Chromatograms of Components Present in Large Quantity in Sample C or Sample D

Fig. 11 shows the mass spectrum of Component 1. The formula prediction results indicated the assignment of  $m/z$  465 as the sodium adduct molecular ion of Bis (triethoxysilylpropyl) monosulfide (TESPM, composition formula:  $C_{18}H_{42}O_6Si_2S$ ), and the assignment of  $m/z$  323.115, 397.190, and 295.084 as fragmentation ions.

Fig. 12 shows the predicted cleavage positions of TESPM. This suggested that cleavage occurred in the same way as in TESPT, shown in Fig. 6.

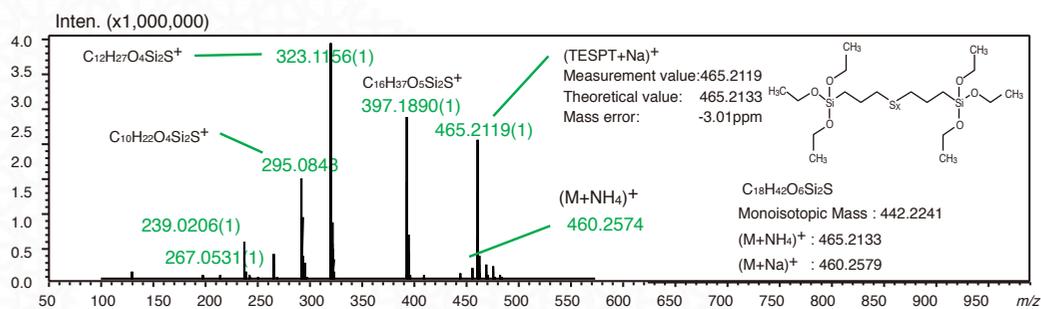


Fig.11 Mass Spectrum of Component 1

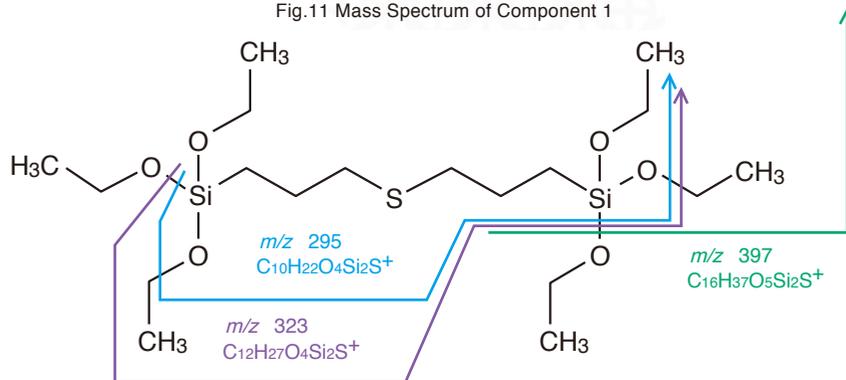


Fig.12 Predicted Cleavage Positions of TESPM

Fig. 13 shows the MS<sup>n</sup> spectra of Component 2. The formula prediction results indicated the assignment of *m/z* 510 as the C<sub>17</sub>H<sub>40</sub>O<sub>6</sub>Si<sub>2</sub>S<sub>3</sub> ammonium adduct molecular ion, suggesting that it had a structure with 1 CH<sub>2</sub> fewer than TESPT. This presents candidates having 2 types of structures, that is, one with the structure where the TESPT -OCH<sub>2</sub>CH<sub>3</sub> replaces -CH<sub>3</sub>, and the other where the methylene chain (CH<sub>2</sub>)<sub>3</sub> replaces (CH<sub>2</sub>)<sub>2</sub>. Next, shifting attention to neutral loss (NL) in MS<sup>2</sup>, where only a neutral loss of 63 was confirmed in TESPT (Fig. 5), in Component 2, a neutral loss of 49 was observed. This corresponds to a decrease of 1 CH<sub>2</sub> in neutral loss as compared to TESPT. Similarly, looking at NL in MS<sup>3</sup>, where only a loss of 74 was confirmed in TESPT (Fig. 5), in Component 2, the neutral loss of 60 was again 1 CH<sub>2</sub> shorter. Thus, the difference is believed to be due to the neutral loss of CH<sub>2</sub> in the structure. As shown in Fig. 6, by comparing the neutral losses it is proposed that the structure of Component 2 is similar to TESPT but in which -CH<sub>2</sub>CH<sub>3</sub> is replaced with -OCH<sub>3</sub>, as shown in Fig. 15.

Formula prediction was also conducted for Component 3, in which it was determined that the formula was C<sub>17</sub>H<sub>40</sub>O<sub>6</sub>Si<sub>2</sub>S<sub>4</sub>. Here we presumed that the structure had 1 S more than Component 2, and even in the sub-components, it was suggested that the distributions are based on the number of S. From the above, Components 2 and 3 are believed to be byproducts or contaminants, confirming that the degree of mixing in such components differed according to the manufacturer.

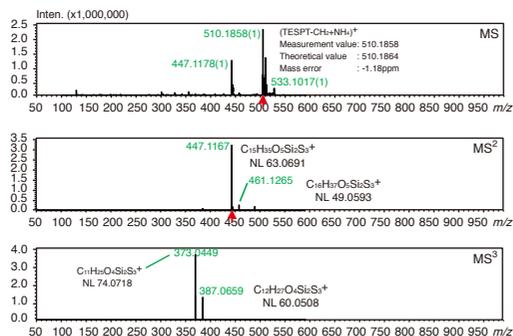


Fig.13 MS<sup>n</sup> Spectra of Component 2

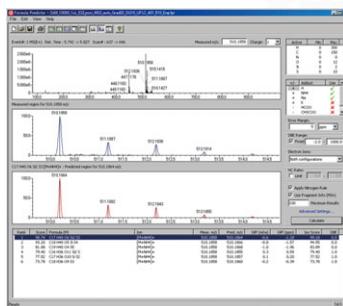


Fig.14 Formula Prediction Results for Component 2

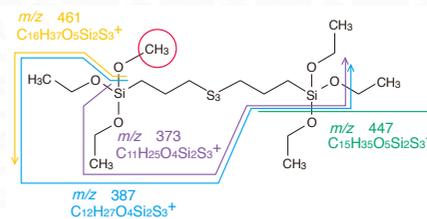


Fig.15 Predicted Structure and Cleavage Positions of Component 2

## 6. Conclusions

- 1) From the TIC, it was observed that the compositions of Sample A and Sample B are similar, and the results supporting this conclusion were obtained through statistical analysis as well.
- 2) It was confirmed from the extracted ion chromatograms that only Sample D was present in TMSPM, and the results supporting this conclusion were obtained through statistical analysis as well.
- 3) Statistical analysis was used to confirm the sub-components that influence Sample C, enabling prediction of their formulas and structures, respectively. This result suggested the possibility that sub-components with distributions based on the number of S, depending on the manufacturer, were mixed in with the products, or the presence of contaminants in the products.

### Acknowledgement

Mr. Fumito Yatsuyanagi and Ms. Yuko Sekine of the Research and Development Center, The Yokohama Rubber Company, Limited gave us the considerable guidance and support for acquiring data.

The information contained in this report is protected by copyright by the publisher, Shimadzu Corporation ("Shimadzu"). The sale, use, reproduction or alteration of this information for any purpose is forbidden without Shimadzu's express written consent, which may be granted or withheld at Shimadzu's sole discretion. Moreover, although the information contained herein has been reviewed, Shimadzu makes no warranty or representation as to its accuracy or completeness. Therefore, if some sort of business or contract is to be entered into in reliance upon this published report and the information contained therein, a Shimadzu business representative must be contacted beforehand for authorization and contractual purposes. This report was created based upon information available to Shimadzu when it was issued, and Shimadzu reserves the right to make revisions to the form and content of this published report without prior notice. First Edition: November, 2010

Founded in 1875, Shimadzu Corporation, a leader in the development of advanced technologies, has a distinguished history of innovation built on the foundation of contributing to society through science and technology. We maintain a global network of sales, service, technical support and applications centers on six continents, and have established long-term relationships with a host of highly trained distributors located in over 100 countries. For information about Shimadzu, and to contact your local office, please visit our Web site at [www.shimadzu.com](http://www.shimadzu.com)



JQA-0376



SHIMADZU CORPORATION. International Marketing Division

3. Kanda-Nishikicho 1-chome, Chiyoda-ku, Tokyo 101-8448, Japan

Phone: 81(3)3219-5641 Fax: 81(3)3219-5710

URL <http://www.shimadzu.com>